
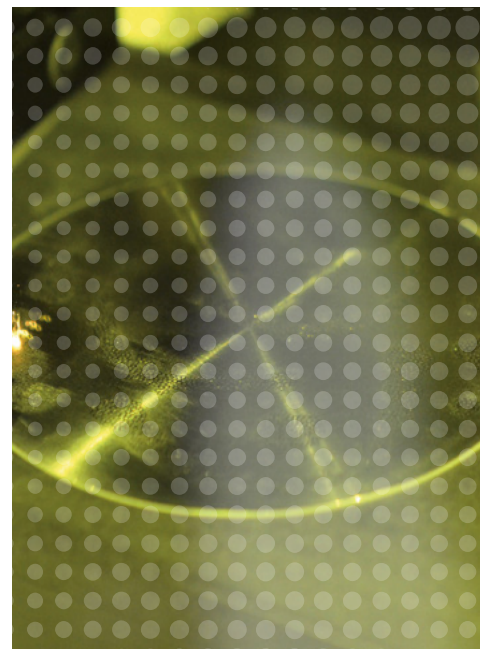



## SMART MEASUREMENT TECHNOLOGY FOR LASERS IN ADDITIVE MANUFACTURING

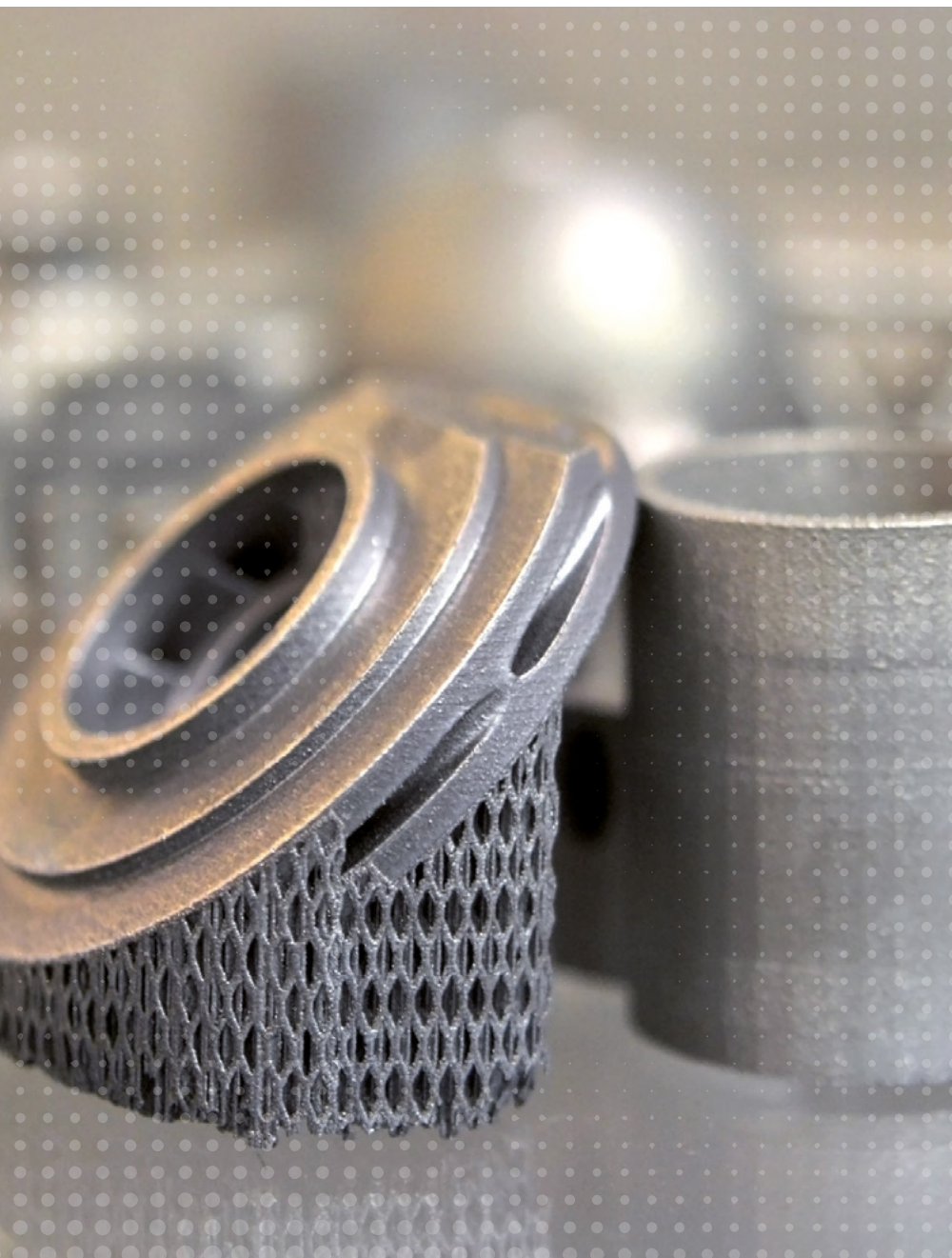

# KNOW YOUR LASER BEAM – ACHIEVE REPRODUCIBILITY


Laser-based additive manufacturing has transformed many industrial processes. But still, one of the most challenging tasks associated with the technology is to ensure reproducibility from one layer to another, between multiple laser sources within the same machine, and from one machine to the next.

The good news is that MKS Instruments has developed a variety of devices especially for additive manufacturing, enabling manufacturers and users alike to measure their laser systems continuously and reliably. We help our customers to record and evaluate ever more laser parameters in the shortest possible time. The Ophir® product portfolio offers the ideal solution for every AM application.

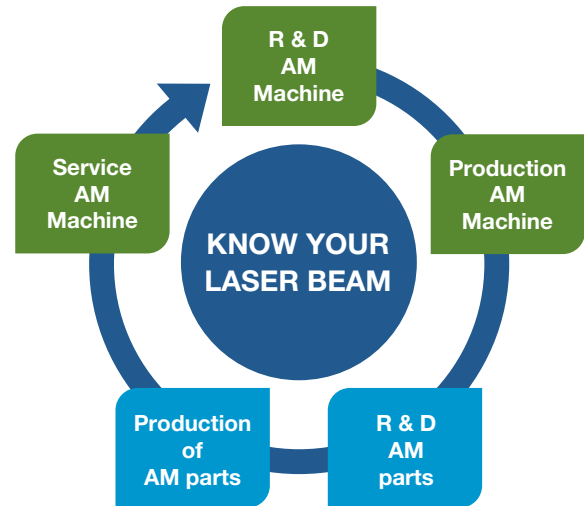

## POWER MEASUREMENT

The first step in knowing whether a laser beam is working according to specifications is to measure the power of the beam. Traditionally, measurements are performed with different types of sensor heads combined with meters or PC interfaces. Designed for use in closed AM production chambers, the compact, robust Ophir Ariel was developed by MKS to work as a standalone power measurement device.

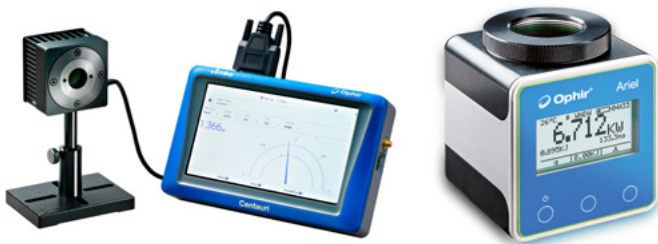


### Laser power measurement

- Determine laser power and energy
- Solutions available for all power levels
- Compact, even with large apertures
- Fast reaction times

## SLIT-BASED BEAM PROFILING

Slit-based beam profiling delivers additional information about the laser beam. The rotational scanning of the slit through the beam creates a 2D image of the beam profile. Our NanoScan beam profilers are compact devices that are ideally suited to measuring low to medium laser power in AM applications. The technology is ISO 11146 compliant.

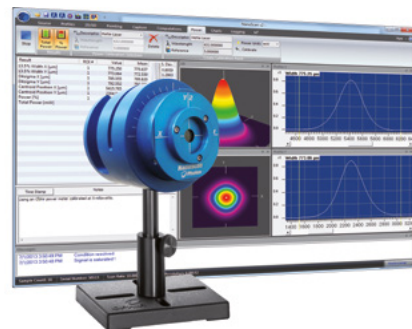


### Slit-based beam profiling

- Determines beam position and size, laser power
- Ideal for additive manufacturing of high-performance plastics
- Extremely compact, suitable for small AM solutions (dentistry, jewelry)

## CAMERA-BASED BEAM PROFILING

Camera-based beam profiling is a popular option for taking a 2D picture of the beam profile. With Ophir BeamGage®, we offer the industry's most advanced beam analysis software; it can be combined with a wide variety of Ophir cameras to fit many applications. Even high-power laser beams can be measured using a CCD camera in conjunction with our Ophir LBS-300-HP-NIR attenuator.

A ready-to-use solution for AM, Ophir BeamCheck combines a high-resolution CCD camera with a calibrated power measurement sensor, allowing for thorough analysis of the intensity distribution within the beam.

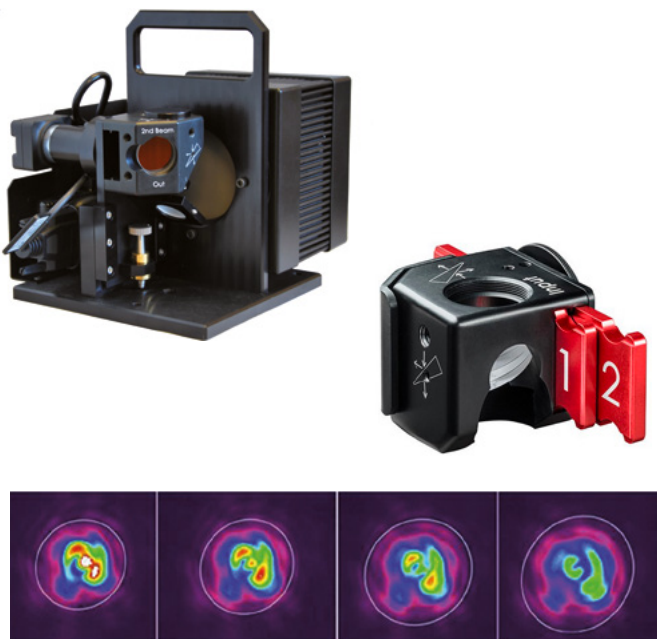


### Camera-based beam profiling

- Focus size, laser power and laser power density
- Short measuring times
- No water cooling required
- BeamCheck additionally includes a calibrated sensor for power measurements

## BEAM CAUSTIC MEASUREMENT

Non-contact beam profiling offers several crucial advantages. First, the measurement technology itself has no power limitation. Even high-power densities upwards of 2 MW/cm<sup>2</sup>, as frequently required in SLM processes, can easily be measured – directly and without attenuation. Secondly, every focus shift is immediately registered and displayed, as a camera detects the Rayleigh scattering of the focused beam in several measurements per second.

In addition, BeamWatch® AM ascertains a variety of beam parameters such as beam position, focus beam size, beam caustics, astigmatism and Rayleigh length. It quickly detects whether the beam is correctly aligned and focused. Powers up to 1kW can be measured over a period of two minutes without the need for active cooling.

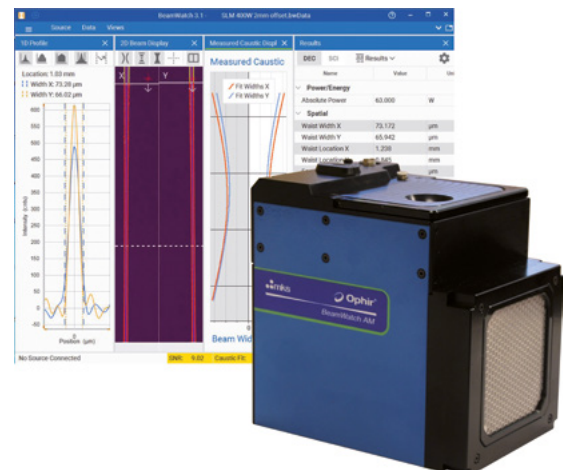

### Beam caustic measurement

- Designed expressly for additive manufacturing
- Non-contact approach to measuring beam profile
- Compact dimensions
- Industrial design with integrated power gauge
- Real-time monitoring of focus shift



## WHY MKS?

<p><b>CRITICAL TECHNOLOGIES</b></p> <p>World-class technology and development capabilities for leading-edge processes</p>	<p><b>PROVEN PARTNER</b></p> <p>Recognized leader delivering innovative, reliable solutions for our customers' most complex problems</p>
	
<p><b>OPERATIONAL EXCELLENCE</b></p> <p>Consistent execution across all aspects of our business</p>	<p><b>COMPREHENSIVE PORTFOLIO</b></p> <p>Largest breadth of product and service solutions for the markets we serve</p>
	

**Ophir-Spiricon, LLC**  
 3050 North 300 West,  
 North Logan UT 84341, USA  
 Tel: +1435-753 3729

**Ophir Spiricon Europe GmbH**  
 Guerickeweg 7, D-64291  
 Darmstadt, Germany  
 Tel: +49-6151-708-0

**Ophir Japan Ltd.**  
 6F Kudan First Place, 4-1-28 Kudan Kita  
 Chiyoda-ku, Tokyo, 102-0073 Japan  
 Tel: +81-3-3556-2781

**Ophir Optronics Solutions Ltd**  
 POB 45021, Har Hotzvim  
 9145001 Jerusalem, Israel  
 Tel: +972-2-548 4444

**WHY OPHIR?** Ophir is a brand within the MKS Instruments Light & Motion division. The Ophir product portfolio consists of high-performance laser and LED measurement technology. Ophir stands for:

- **Stability** – For more than 40 years, Ophir has developed laser measurement systems. This integrates perfectly with the long-term stability and growth of MKS, itself founded in 1961.
- **Variety of products** – Ophir’s product range includes sensors to measure laser power and energy; beam profilers to measure focus shift and beam quality, including industry-leading non-contact measurement systems; and technologies to measure LED luminaires.
- **Individuality** – In addition to the continuously growing portfolio of standard sensors, Ophir develops customer-specific OEM solutions for individual application requirements.
- **Service** – Ophir offers service and calibration centers worldwide that are ISO17025 certified or are in the process of accreditation.

For further information please visit [www.ophiropt.com](http://www.ophiropt.com)